



INSTALLING DRAGON NATURALLY SPEAKING SDK CLIENT 10.1

Pre-Installation Checklist:

1. Refer to the eRAD PACS Workstation Requirements Documentation for the workstation specifications and space requirements to install and run this Speech Recognition Plugin.
2. Remove any previous versions of Nuance Speech Recognition Software. It is recommended to run the Dragon Naturally Speaking Removal Tool if any previous versions of Dragon have been used on the workstation to completely clear any references to the application.
3. Exit all running applications on the workstation.
4. Dragon SDK Client edition requires a DVD drive to install the application.


Installation Instructions:

Place the installation disk in the DVD Drive and the installation should launch. If not, please open the DVD ROM drive in My Computer and open "**setup.exe**"


*Please Note it may take a few moments on this screen to check the system's resources:



Select **Next**:

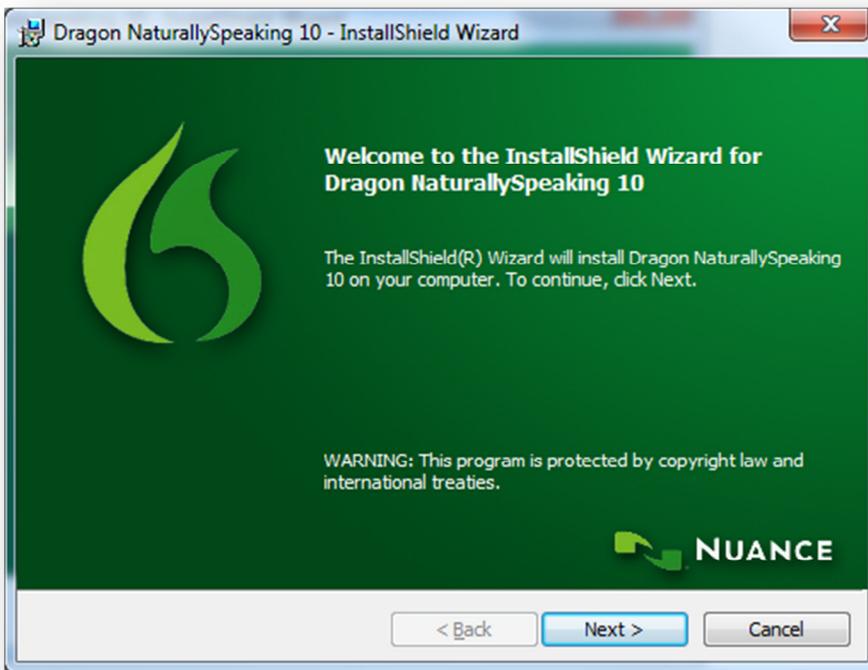

Please Read & Accept the terms in the License Agreement:

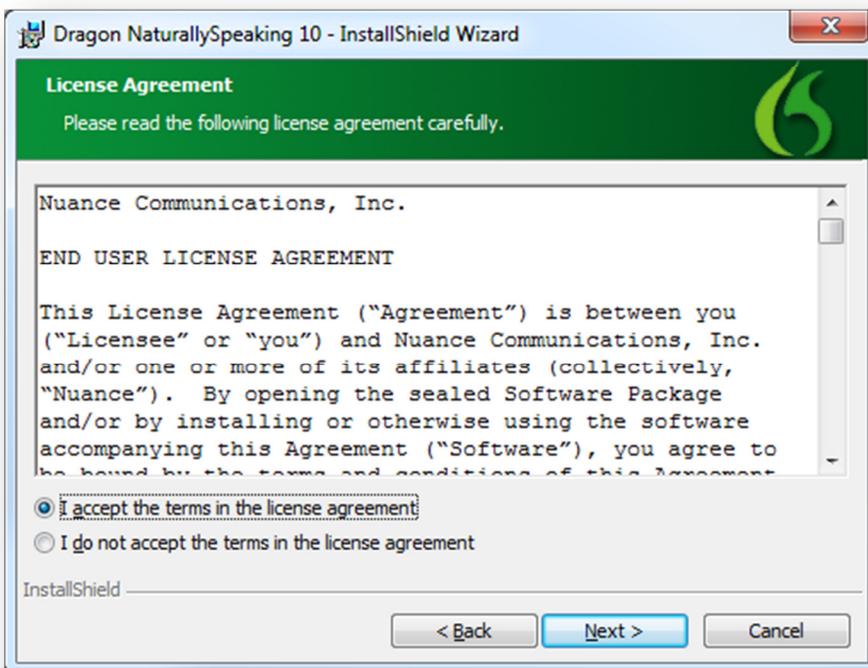

Enter the username and organization, and select **Next**. *The Serial Number will be pre-populated during the installation.

Dragon NaturallySpeaking 10 - InstallShield Wizard

Customer Information
Please enter your information.

User Name: eRAD User


Organization:

Serial Number: - - - - -


InstallShield

< Back Next > Cancel


Select "Custom" installation type and select **Next**. Changing the Installation Directory is optional.

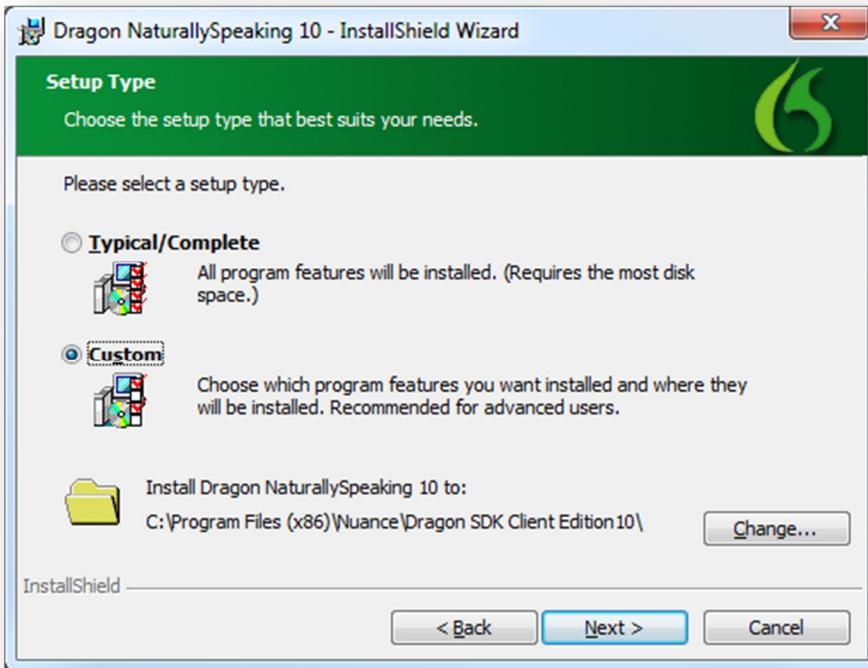


The Custom Setup screen will appear next, ensure you leave it on "run all from my computer" default by clicking the drop-down arrow to the left of "Dragon Naturally Speaking" then select **Next:**

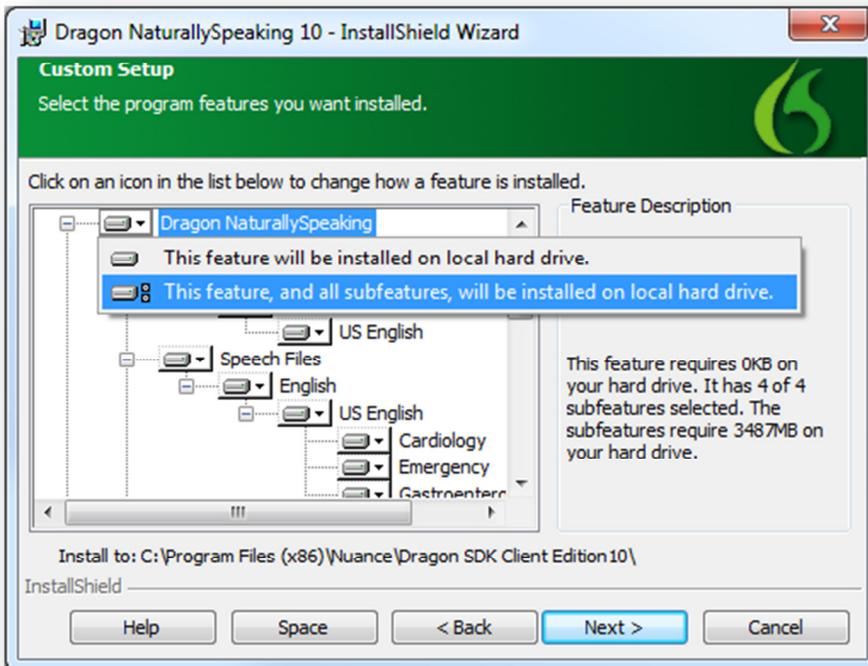


Select the checkboxes to modify the following options during installation and select **Next:**

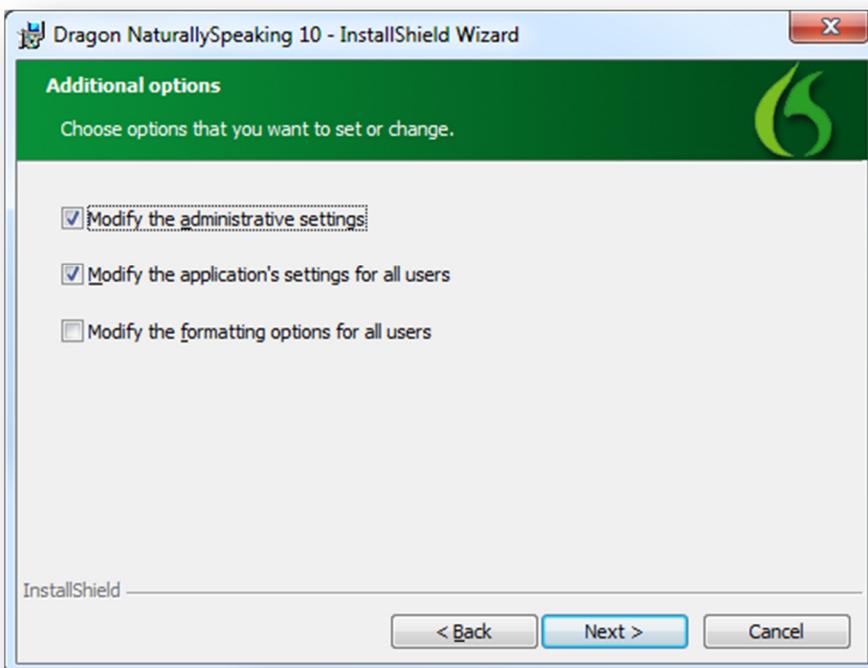


Leave "Enable QuickStart mode" disabled and select **Install**:

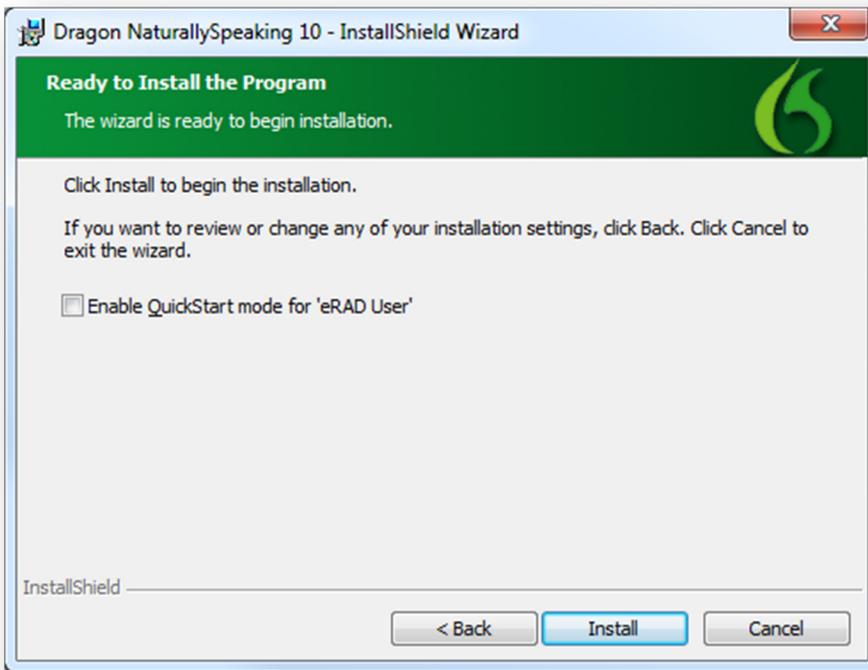

Note: A list of running applications may need to be exited before continuing the installation process by selecting Retry:

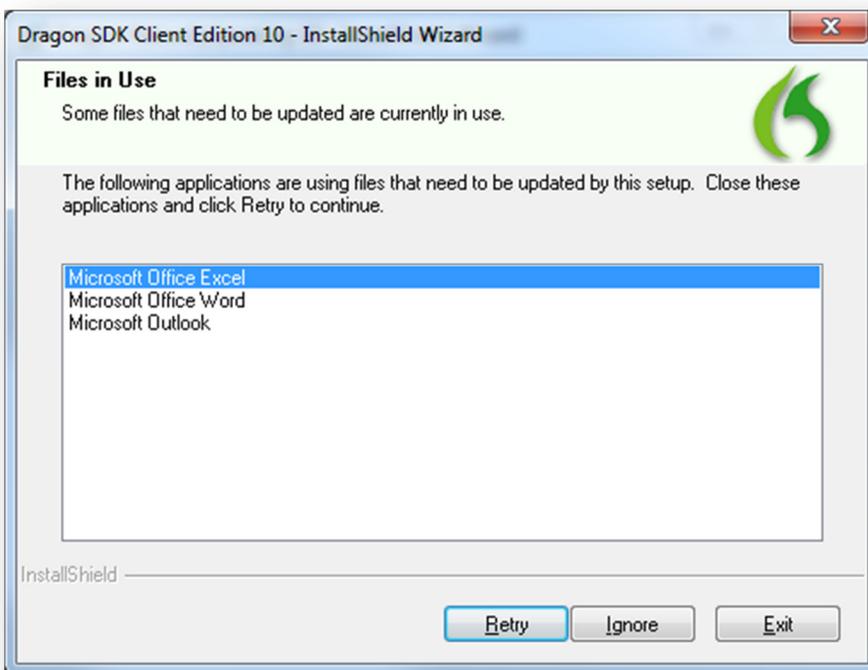


Dragon will now copy files to the workstation. When it is complete it will load the configuration wizard automatically.

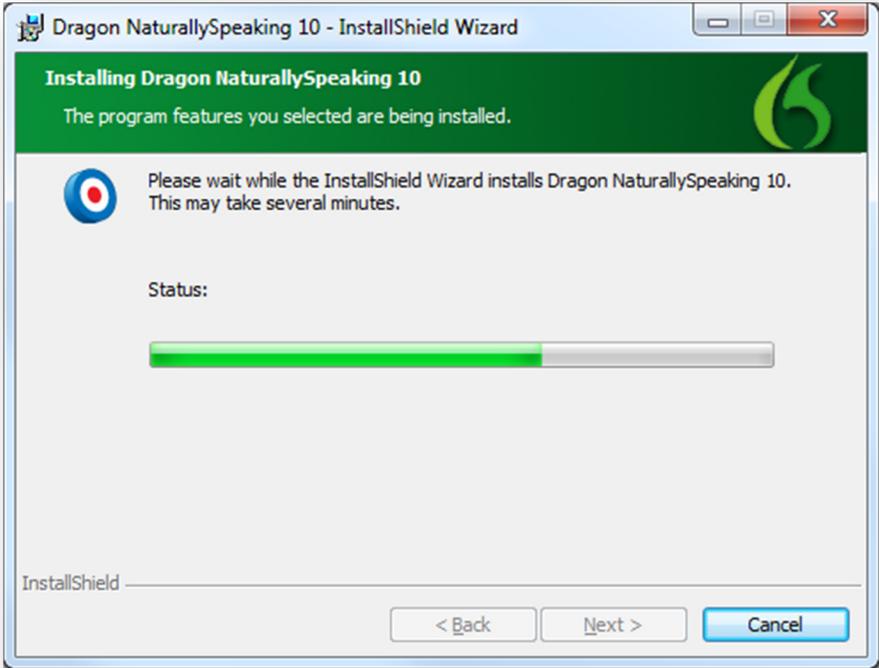


Options Screen Settings

Correction Tab - Confirm only the following options are selected below:

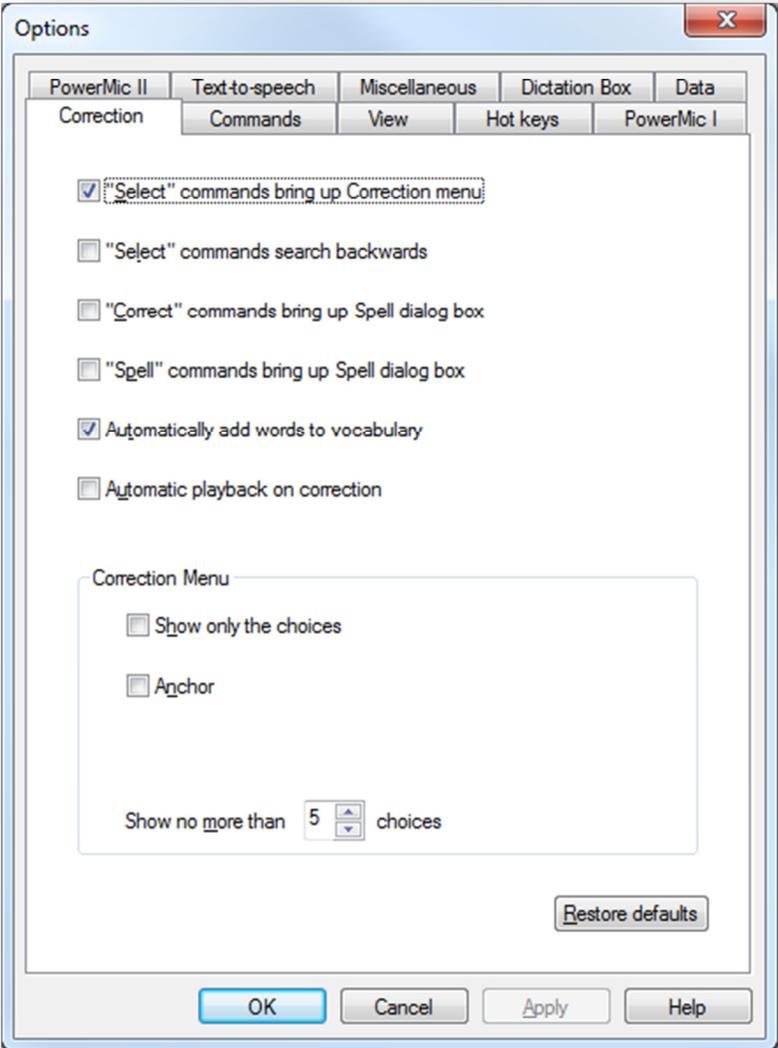


Select **Commands Tab**. Please select the following options:

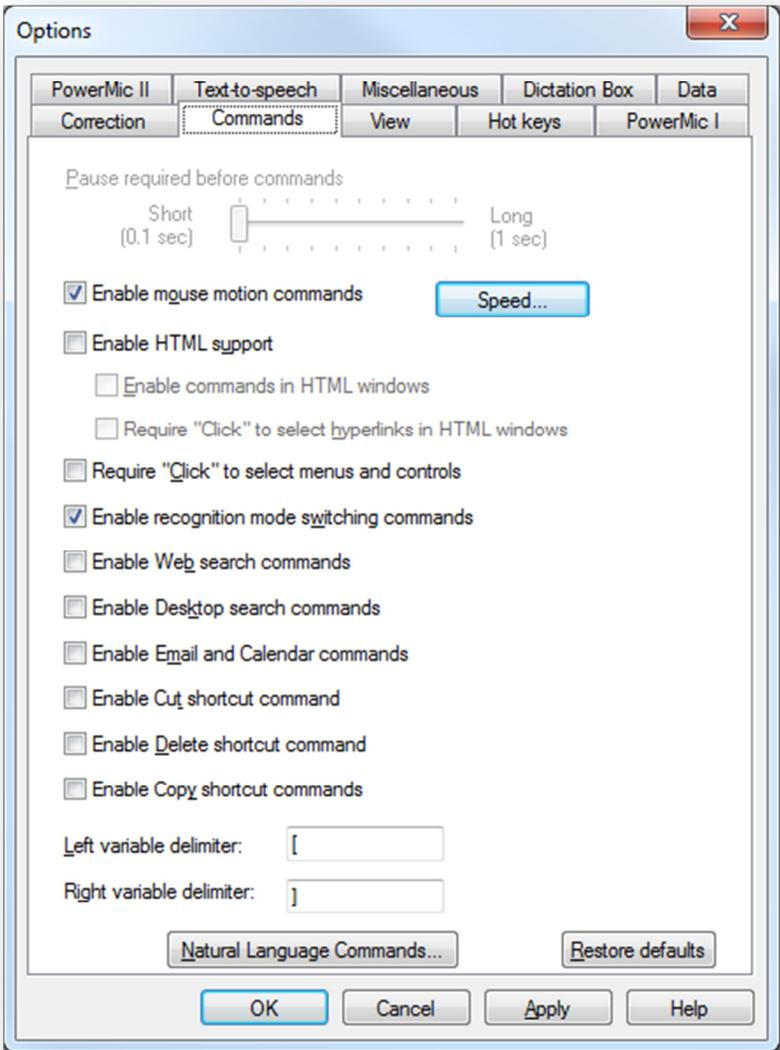


***Note:** These settings may be configured later by individual users to suite their workflow preference.

Select **Data Tab**. Please select the following options, and select **OK**:

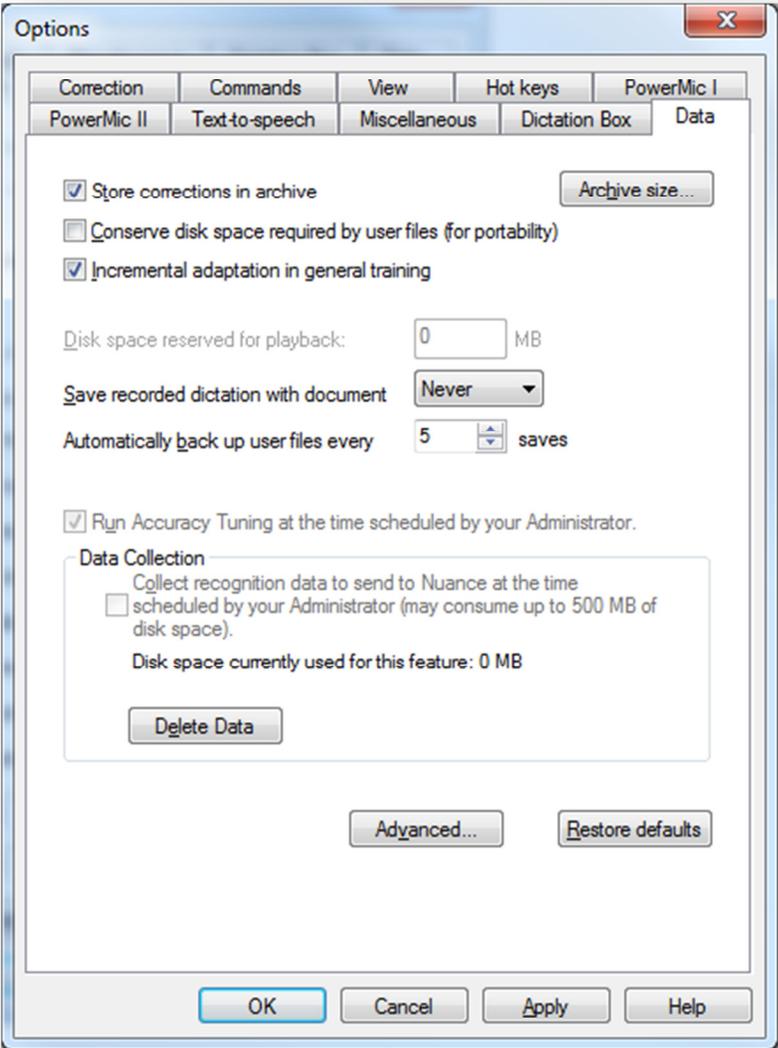


Administrative Settings

Select **Miscellaneous Tab**. Confirm only the following options are selected:

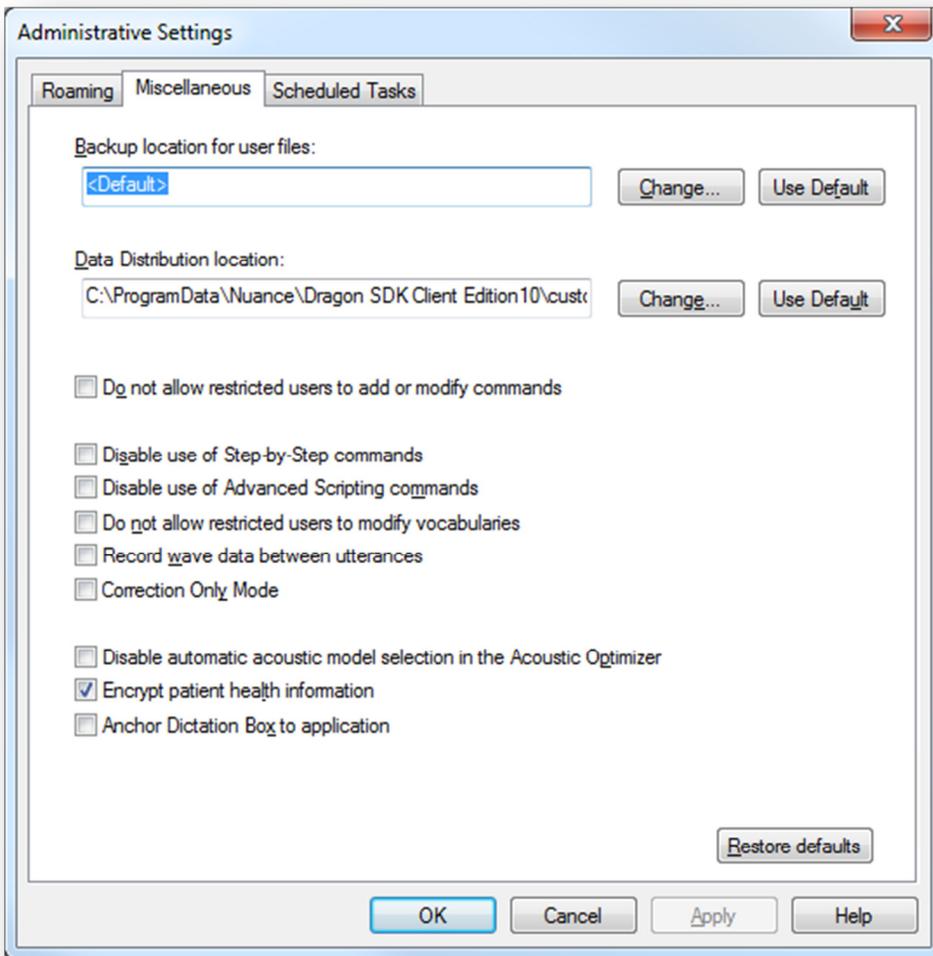

Select **Scheduled Tasks Tab**. Select only the following options and select OK:

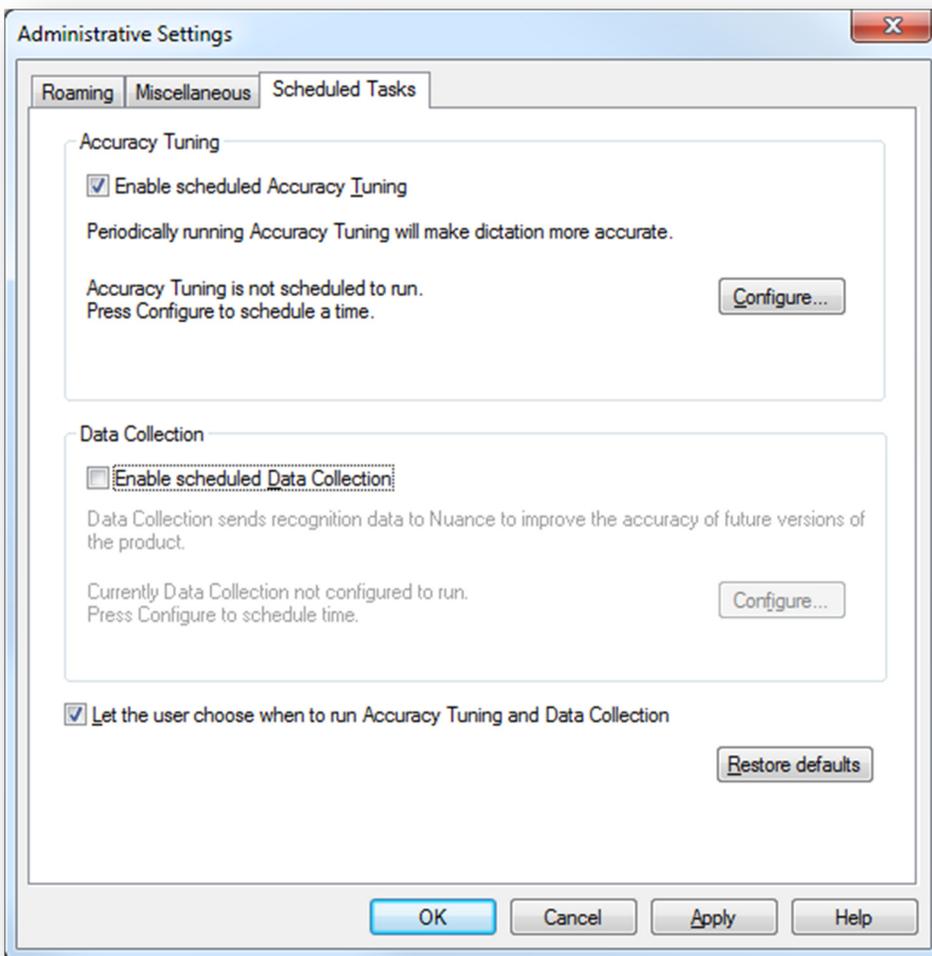

***Note:** To schedule Accuracy Tuning, Windows Administrator rights are required on the workstation login.

Select **Finish**:

